



SUPERCATH 5 & AZ SAFETY IV CATHETERS



THE COMPANY

Global Medics is a family owned company originally established in the 1990's to provide quality product and prompt, professional service to the healthcare market. Ongoing contact with a wide spectrum of medical providers has seen our range expand to include major international brands and our own branded products in four key areas - Aged Care, Theatre/Surgical, Pharmacy/IT and Hospital.

PARTNERS FOR LIFE

Throughout our history, we have strived to forge ongoing relationships with both our customers and our suppliers, allowing us to better meet your needs and provide the best, most efficient and cost effective solutions to improving medical care for your residents, patients and your own healthcare teams.

Truly understanding what you do and working together in this way, we can help enhance quality of life for all patients.

GLOBAL SOLUTIONS Today our range includes a comprehensive supply of medical and healthcare requirements from gloves to wound care, infection control, theatre and surgical items and more that really deliver results.

To view our full product range please visit www.globalmedics.co.nz



OUR MISSION

To forge partners for life so we are better able to deliver the highest quality, innovative products, tested and proven in the global healthcare arena

To reflect the value we place on these partnerships by delivering unbeatable service and backup support

To deliver the most simple and effective solutions that focus on improving quality of life

FOR ALL YOUR MEDICAL PRODUCT ENQUIRIES, BROCHURES OR ORDERS PLEASE CONTACT OUR HELPLINE

0800 GLOMED 45 66 33

Ph 09 920 9900

Fax 09 920 9901

customerservice@globalmedics.co.nz

www.globalmedics.co.nz



Global Medics
Medical and Hygiene Supply Group



INDEX

Global Medics seeks continuously to deliver optimum solutions to your requirements through a range of unique, internationally branded medical products sourced from round the globe.

This brochure outlines our Safety IV Catheter and Dialysis Catheter Range.

COMPANY OVERVIEW	1
INDEX	2
SAFETY IV <i>Needle Stick Protection</i>	3
SAFETY IV <i>Blood Exposure Protection</i>	4
SAFETY IV <i>Notched Needle System</i>	5
SAFETY IV <i>Wings</i>	5
SAFETY IV <i>Features</i>	6
SAFETY DIALYSIS CATHETERS	7
TOURISTRIP <i>Disposable Tourniquets</i>	8
CLINELL <i>IV Wipes</i>	9
<i>Product Codes and Pack Sizes</i>	10




www.globalmedics.co.nz

SAFETY TECHNOLOGY




FULLY ENCASED NEEDLE PROTECTION

- Reduces risk of accidental needlestick injuries to healthcare workers

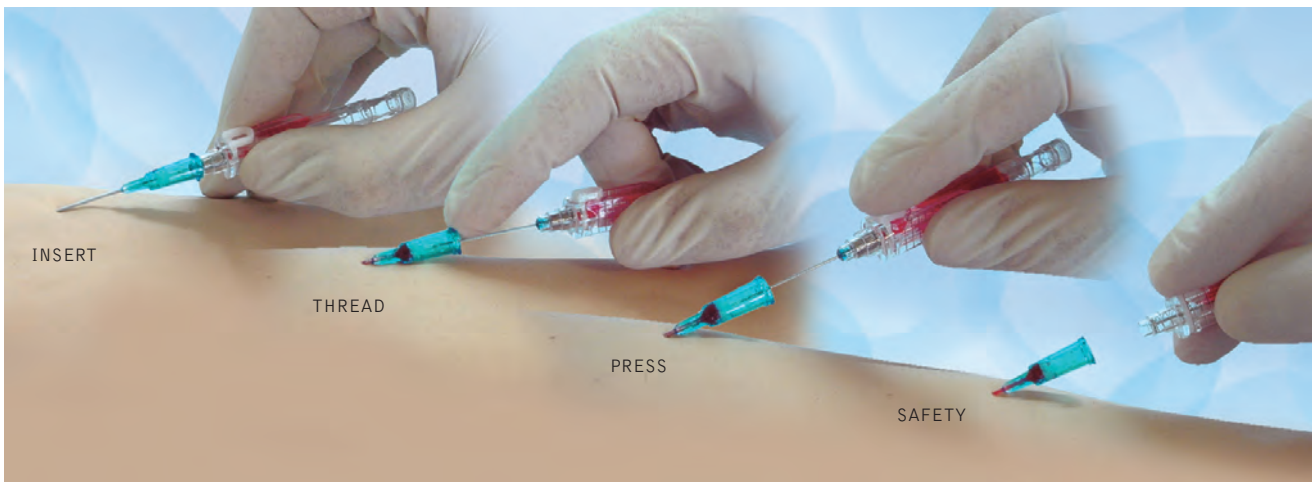


- Protects from blood exposure to contaminated needle
- Offers a reassuring and confident feeling of safety

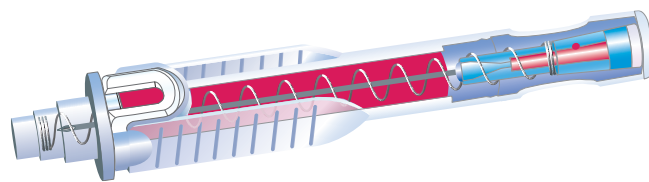

EASY ONE-HANDED SAFETY ACTIVATION

- When the safety mechanism is triggered, introducer needle is quickly retracted into needle shield
- Minimises in-service training


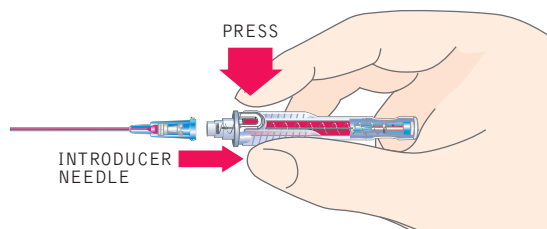


VALVE TECHNOLOGY




PATENTED INTEGRAL CHECK VALVE

- Minimises blood exposure during venipuncture
- No digital pressure required after withdrawal of the needle
- Passive blood exposure protection from the catheter hub

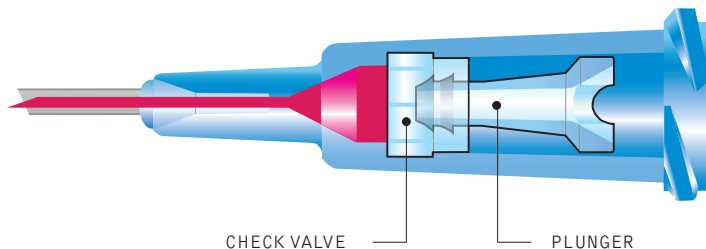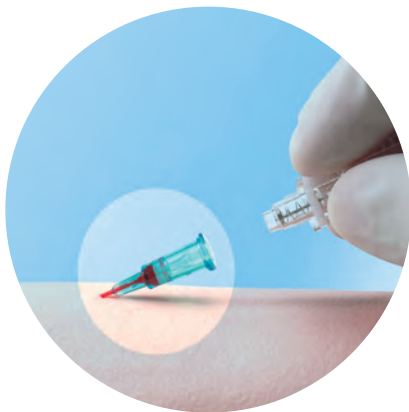

NON-TORTUOUS, STRAIGHT-THROUGH FLUID PATH

- Reduces risk of blood hemolysis
- Ensures a smooth fluid flow

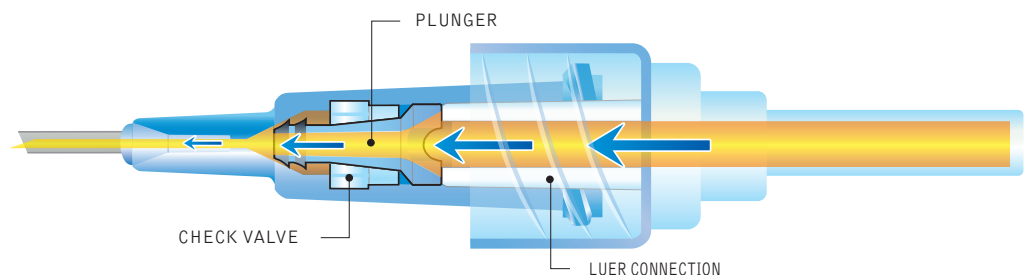

CATHETER HUB DESIGN SAME AS CONVENTIONAL I.V. CATHETERS


- Simpler taping and patient comfort



LATEX FREE & DEHP FREE

- Eliminates dermal reaction to latex and eliminates a substance potentially harmful to patients

NOTCHED NEEDLE SYSTEM FOR QUICK FLASHBACK

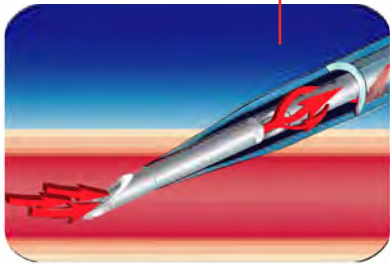


DESIGNED TO ENHANCE SUCCESSFUL VESSEL ENTRY

- Features a side notch in the left side of the needle, just 5 mm to the heel of the bevel

This side notch allows clear and quick visualisation of flashback



- Before Insertion

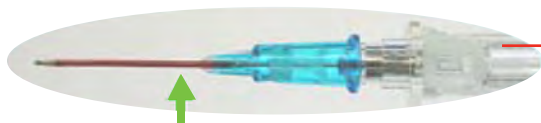

- Immediately After Entry
Flashback from the notch observed between the needle and catheter



- During Insertion Process
Flashback also observed at the flashback chamber



WINGS

An option for those who prefer this style of needle shaft.



PUSH-TO-SHIELD BUTTON

- Easy to activate safety mechanism – indicates needle bevel orientation



ISO STANDARD COLOUR CODING

- Easy identification and embossing of catheter sizes



V-SHAPED NEEDLE POINT

- Ultra-smooth insertion for patient comfort



NOTCHED NEEDLE

- Enhances successful vessel entry



RADIOPAQUE POLYURETHANE

- Enables X-ray detection
- Firm for insertion, softens upon placement for patient safety





CATHETER CONNECTION WITHOUT METAL COMPONENT

- Eliminates potential for false MRI



WINGS

- New option

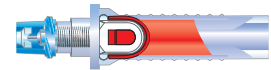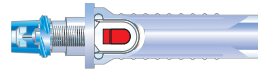

RIBBED GRIP & CLEAR CHAMBER

- Firm grip during procedure
- Fast flashback visibility



NEEDLE SHIELD

- Reduces needlestick injuries
- Blood exposure protection from contaminated needle
- Compact design to save space inside the sharps container



PATENTED INTEGRAL CHECK VALVE

- Minimises blood exposures during venipuncture
- No digital pressure required after withdrawal of the needle from the catheter hub



ROUND CATHETER HUB

- Decreases patient discomfort after securement



FOR DIALYSIS PATIENTS & HEALTHCARE WORKERS

Latest 'Safety' design with fully-encased needle, Supercath AZ helps to prevent needlestick injury and to minimise risk of blood exposure, provides passive safety technology to dialysis patients and healthcare workers.

FEATURES & BENEFITS



A UNIQUE, EASY HANDLING TELESCOPIC NEEDLE RETRACTION SYSTEM

- Eliminates accidental needlestick for both healthcare workers and patients on removal of introducer needle, and prevents blood exposure with fully-encased needle



FEP HD CATHETER

- Minimises risks of internal AV fistulae damage
- Extremely comfortable compared to metal AVF needles
- Highly biocompatible



MALE CAP WITH INTEGRATED STOP VALVE

- No blood exposure after removal of introducer needle



ERGONOMICALLY TAPERED CATHETER

- Catheter tapered up to larger hub diameter enhances flow



LANCET POINT NEEDLE TIP

- Minimises vessel trauma with smooth cannulation and lessens needle patients' discomfort

SINGLE PIECE BODY MOULDING PROCESS

- Reduces hemolysis with maximised flow

Tournistrip was developed by healthcare professionals as a direct result of an unmet clinical need, that being a disposable tourniquet that prevented the cross contamination potential of reusable tourniquets, whilst overcoming the product downsides of other disposable tourniquets.

The result is a unique design, high quality product that is truly single use.







FEATURES & BENEFITS

- Easy to use
- Looks and operates just like a conventional tourniquet
- Easy to tension
- Single handed release and removal

Tournistrip is the only tourniquet to have been clinically trialled. This involved two London teaching hospitals, 226 patients and 20 phlebotomists. The results were overwhelmingly positive:-

- Professional, comfortable and did not pinch their skin
- Phlebotomists saw it as an alternative to their 'gold standard' conventional tourniquets
- 78% preferred Tournistrip to their current tourniquets

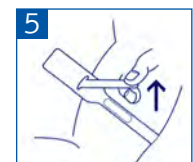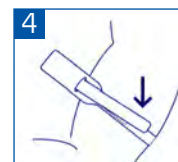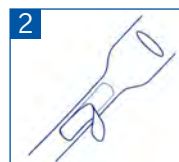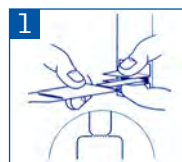

INSTRUCTIONS

1. Remove Tournistrip
2. Peel label to expose adhesive
3. Wrap Tournistrip around arm
4. Tighten and stick Tournistrip down
5. Release Tournistrip





ALCOHOL 2% CHLORHEXIDINE FOR MEDICAL DEVICES & SKIN

Each individual sachet of Clinell Alcoholic 2% Chlorhexidine wipe contains Chlorhexidine Gluconate BP 2% and Isopropyl Alcohol 70% to deliver rapid antimicrobial action and disinfection of medical devices (green pack) and skin (blue pack).



GREEN MEDICAL DEVICE WIPES

- Rapid antimicrobial action with excellent residual activity
- For the disinfection of medical devices such as invasive lines, peripheral intravenous cannulas, renal dialysis catheters and blood culture bottle tops
- The only alcohol 2% Chlorhexidine wipe that is maceratable and biodegradable
- Proven to decrease catheter related bacteraemias
- Only to be used on medical devices



BLUE SKIN WIPES

- Disposable skin disinfection wipes that provide rapid antimicrobial action with excellent residual activity
- For skin disinfection to help reduce bacteria that can cause skin infections
- The only alcoholic 2% chlorhexidine wipe for skin on the market

STUDY:

USE OF CLINELL ALCOHOLIC 2% CHLORHEXIDINE ON CENTRAL VENOUS CATHETER PORTS

This was done at Great Ormond Street by Dr James Soothill, Microbiologist, with the aim of seeing if using these wipes would reduce the incidence of infections through central lines. It was implemented on 5 wards and used for 6 months. The rate of line infections was compared to the previous 6 months in terms of infection rates.

RESULT:








- All 5 wards showed a drop in line infection rates
- All 5 wards showed a decrease in positive blood cultures





PRODUCT CODES

SAFETY DIALYSIS CATHETERS **Supercath™ 5**

COLOUR CODE ISO10555-5:1996	CATALOGUE NUMBER	GAUGE	CATHETER					PACKAGING	
			FLOW RATE mL/min	O.D	I.D	LENGTH		PCS/BOX	BOXS per carton
 ORANGE *	SP125-14-31	14G	272	2.1 (mm)	1.60 (mm)	1 ¼" (31mm)	NO CHECK VALVE	50	20
 MEDIUM GREY *	SP125-16-31	16G	183	1.7 (mm)	1.25 (mm)	1 ¼" (31mm)	NO CHECK VALVE	50	20
 DEEP GREEN	SP120-18-31	18G	100	1.3 (mm)	0.9 (mm)	1 ¼" (31mm)		50	20
 PINK	SP120-20-31	20G	55	1.1 (mm)	0.7 (mm)	1 ¼" (31mm)		50	20
 DEEP BLUE	SP120-22-25	22G	31	0.9 (mm)	0.6 (mm)	1" (25mm)	OPTIONAL NOTCHED NEEDLE	50	20
 DEEP BLUE	SP120-22-31	22G	28	0.9 (mm)	0.6 (mm)	1 ¼" (31mm)	OPTIONAL NOTCHED NEEDLE	50	20
 YELLOW	SP120-24-19	24G	16	0.7 (mm)	0.4 (mm)	¾" (19mm)	OPTIONAL NOTCHED NEEDLE	50	20

Supercath™AZ

COLOUR CODE ISO10555-5:1996	CATALOGUE NUMBER	INNER NEEDLE GAUGE	CATHETER					PACKAGING	
			FLOW RATE mL/min	O.D	I.D	LENGTH		PCS/BOX	BOXS per carton
GREEN	SP302-16(25)	16G	N/A	2.10 (14G)	1.60 (16G)	25		100	10
GREEN	SP302-16(33)	16G	N/A	2.10 (14G)	1.60 (16G)	33		100	10
ORANGE	SP302-17(25)	17G	N/A	1.93 (15G)	1.45 (17G)	25		100	10
ORANGE	SP302-17(33)	17G	N/A	1.93 (15G)	1.45 (17G)	33		100	10
PINK	SP302-18(25)	18G	N/A	1.73 (16G)	1.25 (18G)	25		100	10
PINK	SP302-18(33)	18G	N/A	1.73 (16G)	1.25 (18G)	33		100	10

TOURNISTRIP DISPOSABLE TOURNIQUETS



COLOUR	CATALOGUE NUMBER				PACKAGING	
					PCS/BOX	BOXS per carton
WHITE	TOURNI1				200	1

CHLORHEXIDINE WIPES



COLOUR CODE	CATALOGUE NUMBER	WIPE SIZE	PACKAGING	
			PCS/BOX	BOXS per carton
GREEN	CA2C20 - 0 Chlorhexidene IV Wipes	19 X 10.5 cm	200 sachets	24 x 200
BLUE	CA2CSKIN Chlorhexidene Skin Wipes	10 X 10 cm	200 sachets	24 x 200

* Please refer to the package insert for complete instructions for use



SUPERCATH 5 & AZ SAFETY IV CATHETERS

FREEPHONE 0800 GLOMED 45 66 33
PHONE 09 920 9900 | FAX 09 920 9901
EMAIL customerservice@globalmedics.co.nz
WEB www.globalmedics.co.nz

